


iM3 Vital Signs Monitor



iM3
8" Color TFT



Stunning Design

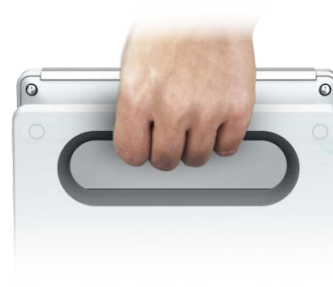

Comfortable handle



Ultra Slim

<2.5kg



EDAN
A world of potential

Full-touch Operation

An 8 inch **capacitive touch screen** is applied in iM3 vital signs monitor which brings the visual enjoyment and highly responsive operation experience.



8 inch
Touch Screen



One button operation

Besides full-touch screen, the waterproof silicon buttons and control knob are integrated in iM3 for a rapidly operation.



■ NIBP

Press the button to start or stop NIBP measurement.

■ New patient

Press the button to admin or add a new patient.

■ Patient Type

Press the button to switch the patient type.

■ Menu

Press the button to enter main menu or back to main screen.




■ Control knob

Easy to use




Rich Display Information

Different from traditional vital signs monitor, iM3 comes with a large screen which enable rich information display. Shortcuts and arrow key are provided for switching among standard display, trend display , trend review and alarm review.

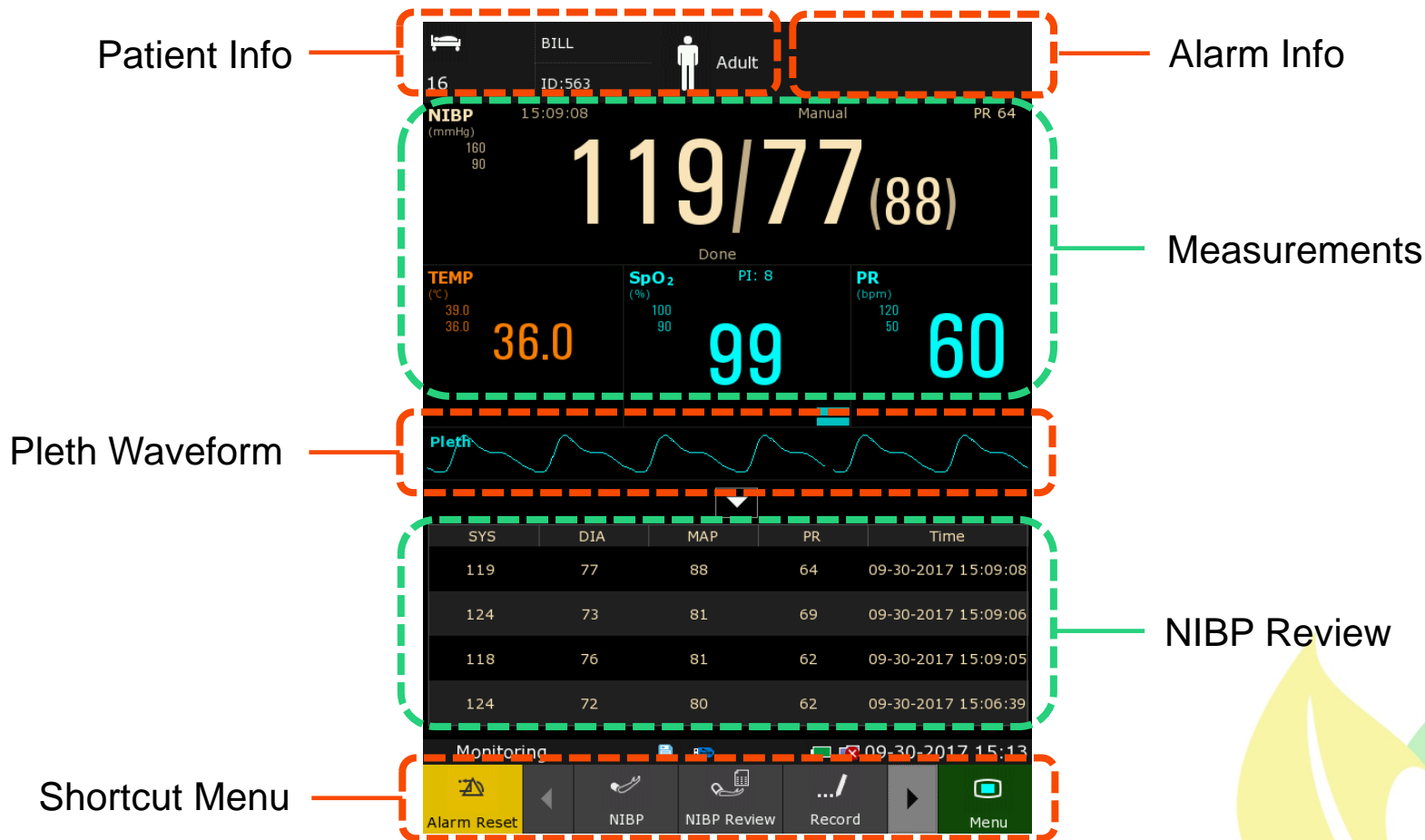

Customizable Shortcut Menu

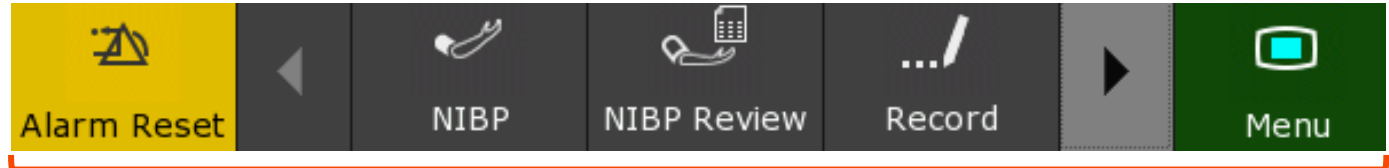

The shortcut menu is put on the bottom of the display, providing direct access to daily operations. And the shortcuts are shifty according to the work mode.

■ Monitoring Mode

Alarm Reset / **Alarm Mute** / NIBP / **Review** (Trend Graph, Trend Table, Alarm, NIBP, History Patient) / Record / General / Standby / Night Mode / Score / Barcode / Admission / Menu

■ Spot Mode

New Patient / NIBP / Barcode / Record / **NIBP AVG** / Review / Standby / General / Score / Menu

■ Round Mode

Select Patient / NIBP / Barcode / Review / **NIBP AVG** / Standby / General / Score / Menu

